

## A generic Microfluidic Approach for Deciphering Nanoscale biovesiclES propertieS

Lyne Pillemont, Daniel Guneysu, Céline Elie-Caille, Wilfrid Boireau, Anne

Marie Gué

### ► To cite this version:

Lyne Pillemont, Daniel Guneysu, Céline Elie-Caille, Wilfrid Boireau, Anne Marie Gué. A generic Microfluidic Approach for Deciphering Nanoscale biovesiclES propertieS. 3rd Caparica Christmas Conference on sample Treatment, Dec 2018, Caparica, Portugal. 2018. hal-02473108

## HAL Id: hal-02473108 https://laas.hal.science/hal-02473108

Submitted on 10 Feb 2020  $\,$ 

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers. L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



# A generic Microfluidic Approach for Deciphering Nanoscale biovesiclES propertieS

PILLEMONT Lyne<sup>1,</sup>, GUNEYSU Daniel<sup>2</sup>, ELIE-CAILLE Céline<sup>2</sup>, BOIREAU Wilfrid<sup>2</sup>, GUE Anne-Marie<sup>1</sup>

lpillemo@laas.fr


<sup>1</sup> LAAS-CNRS, 7 avenue du colonel Roche, F-31400 Toulouse, France <sup>2</sup> FEMTO-ST, 15B avenue des Montboucons, F- 25030 Besançon, France

INTRODUCTION

Extracellular vesicles (EVs) are subcellular species present in all biological fluids which are increasingly considered as pathologies indicators. Therefore identifying and characterizing nanoscale vesicles in biofluids could lead to unprecedented diagnosis/prognosis assays for a wide set of pathologies. Owe to their complexity in size, origin, membrane markers, there is currently no ideal technology available yet to quantify and characterize cell-derived microvesicles, and fully relate their structure and functions which have high potential in clinical diagnostics and early treatment of life threatening diseases such as cancer, cardio vascular diseases....

In that goal, we propose to couple a microfluidic device (lab-on chip) enabling to separate particles from a complex sample to a second device allowing the capture of particles on specific spots and their analysis. These methods allow obtaining important information about nanoscale vesicles in terms of size, membrane markers expression, and structure-function properties.

<b><u>1- Separation Principle</u>: Hydrodynamic filtration</b>	<u>2- Device Design</u> :			<b><u>3-Microfabrication</u>: Lab on chip in clean room</b>
Exclusion of large particles from the streamlines located close to a wall	r <sub>c</sub> and microchannel geometry (width, height) defined:			Fabrication of multilevel microfluidic devices by combining a
Critical continue $\pi_{-}$ - portion of fluid Oin Oout	$r_c = 350$ nm	ain channel Output	Lateral channel	lithography technique and lamination technology

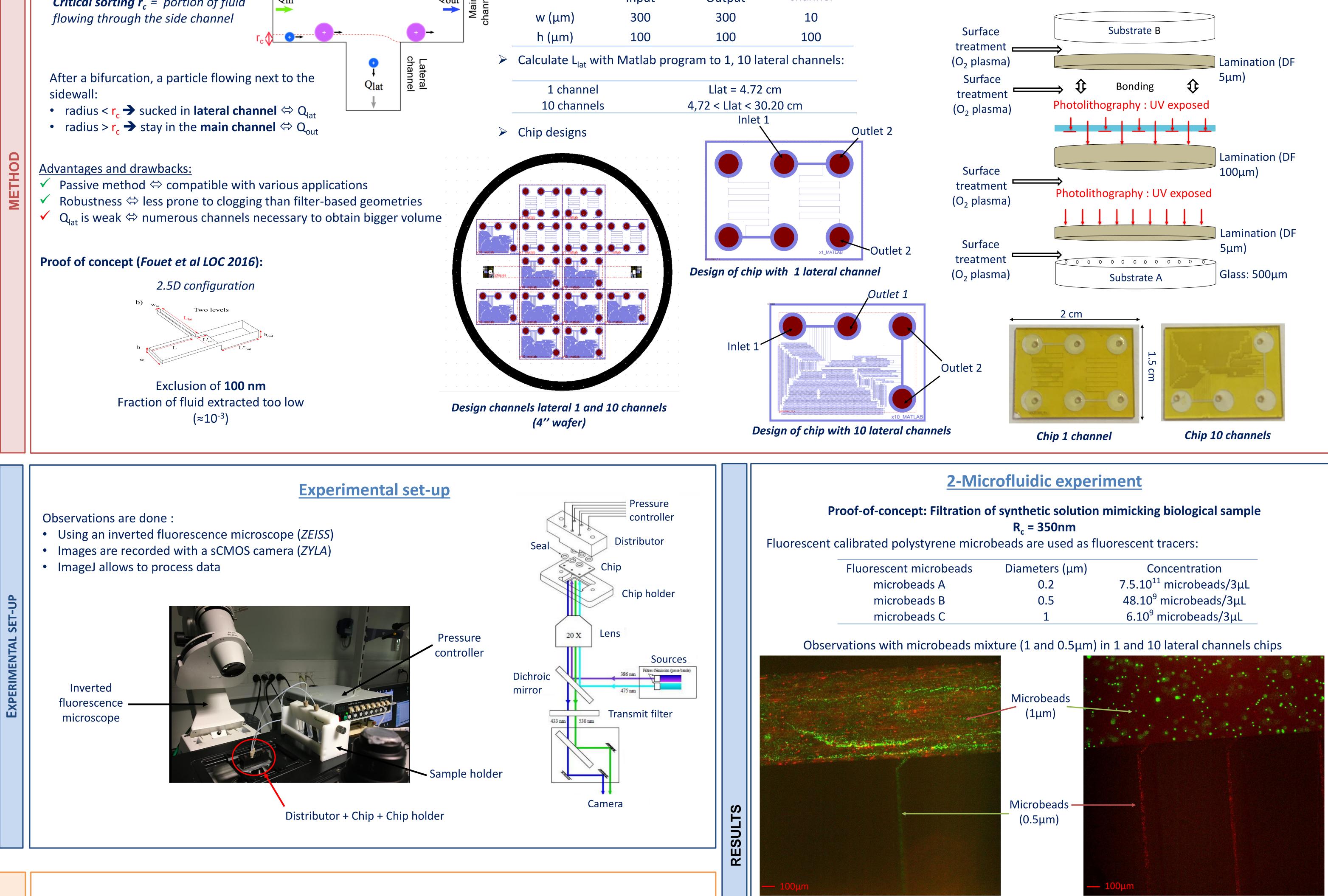






Fluorescent microbeads	Diameters (µm)	Concentration
microbeads A	0.2	7.5.10 <sup>11</sup> microbeads/3µL
microbeads B	0.5	48.10 <sup>9</sup> microbeads/3µL
microbeads C	1	6.10 <sup>9</sup> microbeads/3μL

The critical radius is checked in different chips: the 1µm microbeads cannot flow into lateral channels and 0.5µm microbeads through in lateral channels

# Observations with microbeads mixture (0.2 and 0.5µm) in 10 lateral channels chip


In 10 lateral channels chips, 0.2 and 0.5µm microbeads flow into

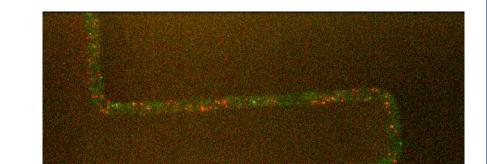

lateral channels

Filtration of synthetic solution on different chips	Mi
confirm the correct critical radius r <sub>c</sub> =350nm	(
c	NЛi

Channels	Flow rate (µL/s)
1	0,000014
10	0,023
10	0,023

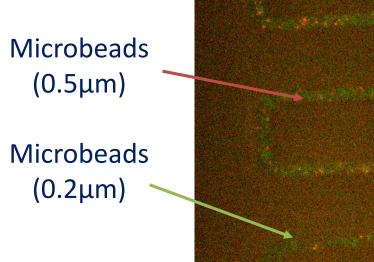

The purpose of this project is to develop a comprehensive microfluidic method which allows micro and nanoparticles sorting by size (hydrodynamic filtration) and to combine this device with a second device which allows capturing these particles on specific spots.

The proof of concept has been done with synthetic solution using fluorescent calibrated polystyrene microbeads.

# **PERSPECTIVES:**

CONCLUSION

- Optimization of filtered volumes (with several viscoelastic effects)
- Increase number of lateral channels
- Filtration test with blood sample
- Make new design with critical radius smaller for exosome applications

	1	0,00014				
	10	0,023	30µm			
ACKNOWLEDGEMENTS						
ANR	Madness projet					
Fren	ch RENATECH network					
REFE	RENCES					
M.Fouet, MA. Mader, S.Iraïn, Z.Yanha, A.Naillon, S.Cargou, AM. Gué and P. Joseph, Lab chip, 2016						







